Mixer Guide

Caframo. LabSolutions





654 Petrolia Road, Toronto, ON M3J 2W3 Tel: (416) 736-6166 Toll-Free: (800) 387 3570 Fax: (416) 736-9346

The Crossover

- Make your job easy A fully featured mixer Powerful, portable plug-in installation
- Capture data Speed, torque, time, date
- Mix batches Media, buffers, creams
 No motor burn out
- Build into a modular system and control by a PLC - LabVIEW compatible
- Safety certified to applicable standards

Built for Real Work

- High efficiency, 300 watts at shaft
- Brushless DC motor fan cooled
- Overload and overheat protected



	1540		
Speed rpm	50-1500		
Volume L (gal)	200 (55)		
Viscosity cP	180,000		
Torque Ncm	3,000		
Horsepower	1/2		
Weight kg (lbs) 12 (26.5)			
Voltage	120 volt		
Work Smart	Digital Data Collection		

O MAX 2500

00:00:00

Includes: 1540 mixer, chuck, key and guard, screen protector, USB, mounting hardware and tools, grounded cord, manual, calibration certificate.

Auto Start AutomationThrough Shaft Adjustment

• Timer Notifications

• Rotation Selection

Three year warranty.

"A powerful and capable stirrer...
Our Techs are quite impressed"


Selection Menu

Set-up Accessories

GUI Resistive Touch Control

- Operate with a gloved hand
- Large visibility of settings
- Speed control on home screen

Settings Menu

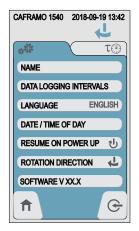


Torque Timer Menu



For full specifications

CaframoLabSolutions.com

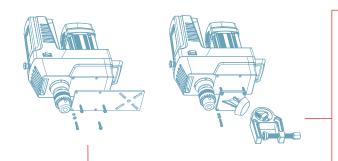

A129 Plate Mount

- 260 x 100 mm (10 x 4 in)
- 6.4 mm (¼ in) thickness
- 304 electropolished stainless steel

Sturdy simple mounting solution. Secure to existing horizontal surfaces, mixing equipment, drum brackets or industrial stands. Slotted bolt pattern provides adjustability when mounting.



Download template for mounting.

A128SET C-Clamp Set

- Single bolt adjusts angle on two axes
- Attaches directly to a drum, frame or stand
- Cast 304 electropolished stainless steel

Mount the Caframo Crossover mixer to any stable vertical including an industrial stand, drum or vessel. Easy single bolt adjustment for multi-angle positioning.
Set includes c-clamp, clamp adapter and

Set includes c-clamp, clamp adapter and hardware. Standard tool installation.

A118 Stand with Casters

- Adjustable 99 to 145 cm (39 to 57 in) height
- 102 x 102 cm (40 x 40 in) caster footprint
- Plate steel, square tube construction


This manual lift, heavy duty stand and base provides a solid mount for the Caframo Crossover when mixing in drums or tanks up to 122 cm (48 in) in diameter. This stand comes with smooth rolling casters that swivel and lock in place.



Accessories



Communication Ports

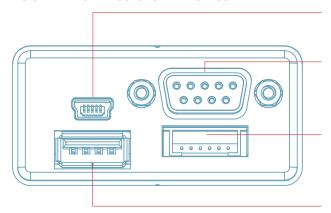

Mini USB

PC Control - communicate using command set.

DB9 Serial Port - RS232

PC or PLC Control and legacy ethernet connectivity. Adapt to RS485 and daisy chain mixers together.

Analog 4-20 mA

Stable PLC control over long distances.
Simple wiring scheme widely used in industry.

Standard USB

Data logging and software updates. Process starts as soon as memory stick is inserted.







CG6	A354	A356	A016SET
Chuck Shaft Guard	Propeller with Shaft	Propeller with Shaft	Chuck Set
Polycarbonate ABSSimple mountingSafety accessory	 102 (4) dia blade 914 x 12.5 (36 x 1/2) shaft 316L SS welded and polished 	 152 (6) dia blade 914 x 12.5 (36 x 1/2) shaft 316L SS welded and polished 	Shafts up to 16 (5/8) diaThrough shaft420 SS

mm (in)

800-567-3556 ext 248
CaframoLabSolutions.com
501273 Grey Road 1, Georgian Bluffs ON N0H 2T0

